

MZ-dent[®]
GERMANY

Low Speed Handpieces



Discover the advantages...

Discover the advantages

The MK-dent straight and contra angle handpieces offer high tech precision and reliability for attractive prices. First have a look at the contra-angle handpieces. You can choose between LS2011L with glass rod light and the non optic version LS2011. Take it in your hand and feel the excellent grip.




25.000 LUX glass rod



Full compatible to all ISO-standard/INTRAMatic® motors



Good grip, relaxed and safe work at all times

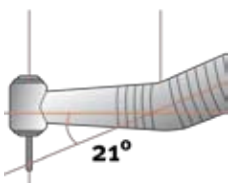

Special contra-angle design prevents tilting of the instrument

For good reasons the instrument is always very precise to handle. It has the right weight, not a gram too heavy or a gram too light. During preparation, the pressure point on the tooth is exactly on the extrapolated axis of the contra-angle handpiece and thus prevents tilting of the instrument. The two-part design makes it economic and flexible. Handpiece and the new MK-dent instrument head LS9011 and LS9021 are quick, easy and safe connected.

Further features of the heads LS9011/LS9021:



- 1:1 transmission LS9011 / 2:1 transmission LS9021
- for removing carious material
- maximum operating speed 40.000 rpm LS9011 / 20.000 rpm LS9021
- dirt guard and strengthened chuck system
- for changing of instruments you just push a button

And as easy is the connection of the contra-angle handpieces to the supply system just by a Click. They are compatible to all ISO-standard/INTRAMatic®. Comparable benefits you have with the new Straight Handpiece LS1011 from MK-dent. Like the contra-angle types with separate air and water supply lines incorporated. Low noise and vibration during operation. High quality stainless steel, easily sterilizable in an autoclave at 135° C, easily washable in a thermal disinfectant. Very long lifetime.



AM1014 Airmotor

- up to 20.000 rpm
- infinite right and left handed rotation adjustment
- constant speed
- separate cooling mediums
- long lifespan because of high quality spare parts
- easy to clean and care
- fits for MK-dent LS2011 and LS1011
- KaVo* compatible

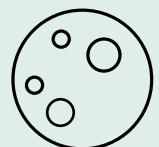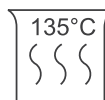



- ① **LS2011L** contra angle handpiece with light
Body 1:1 transmission, blue, 25.000 LUX fibre optic glass rod, max. 40.000 rpm, compatible to all ISO-standard/INTRAmatic® e- and airmotors, KaVo* compatible
- ② **LS2011** contra angle handpiece without light
Body 1:1 transmission, blue, max. 40.000 rpm, compatible to all ISO-standard/INTRAmatic® e- and airmotors, KaVo* compatible
- ③ **LS1011** straight handpiece without Light
Body 1:1 transmission, blue, max. 40.000 rpm, compatible to all ISO-standard/INTRAmatic® e- and airmotors, KaVo* compatible
- ④ **LS9011** contra angle head
Head 1:1 transmission, blue, max. 40.000 rpm, fits to MK-dent LS2011L, LS2011 and KaVo* contra angle handpieces, comfort pushbutton system
- ⑤ **LS9021** contra angle head
Head 2:1, green, max. 20.000 rpm, fits to KaVo* contra angle handpieces 2:1 and 7,4:1 transmission, comfort pushbutton system

MK-dent[®]

G E R M A N Y



MK-dent is committed to sell high quality dental equipment and supplies through authorized dealer network 's worldwide. More about us and our program, including Low Speed Handpieces, spare parts, tools, care products and more, you'll find at our website: www.mk-dent.com, in our special catalogues, data sheets and of course by personal information through our dealers and German headquarter.

Authorized Dealer: